




TE II Lab Centers

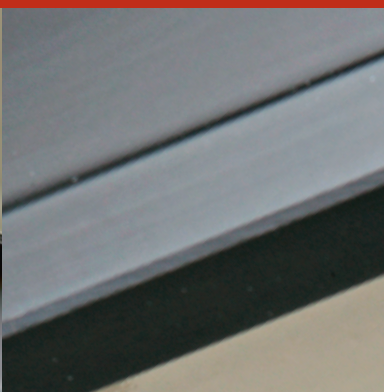

TE II Lab Centers

Total Experience® II Lab Centers integrate the work counter, standing area and traffic lanes, using every inch of available floor space to maximum benefit. Sheldon's unique TE II design eliminates the need for students to work back-to-back, creating a spacious, safe working environment.



The TE II allows each student ample front-counter space and elbowroom. Other workers and activities are always within the student's angle of vision, making lab work not only easier but safer. Each student can work with either neighbor across from his or her corner, and because students are not working back-to-back and obstructing aisles, traffic lanes are clear in every direction.

Available in a variety of surface finishes, the TE II is ideal for chemistry and physics labs and is well-suited for biology and general science applications.

"The TE II Lab Centers free teachers from the standard peninsular-style laboratory design. These centers assist with flow and safety for students and teachers."

James W. Collins
Science Program
Coordinator
Charles A. Dana Center
The University of Texas
at Austin

Features

The TE II is available in three surface materials:

The **epoxy-resin** surface is the most chemically resistant surface available and is ideally suited for chemistry applications.

The **electronically laminated, edgegrain-maple** surface provides an excellent work surface for physics applications.

The **phenolic-resin** surface is chemically resistant and an excellent choice for chemistry, physics, biology and multi-science applications.

The understructure on all models is fashioned of sturdy, molded fiberglass bolted to a 10-gauge steel support pedestal. Two removable panels provide easy access to plumbing.


Centers are equipped with two duplex electrical receptacles and are designed to accommodate optional storage cubicles.

The TE II is available in 17 models featuring varying surface materials and tabletop sizes. Some models are equipped with round or trough sinks and Sheldon's vandal-resistant Unicast fixtures.

Seven TE II models are ADA Accessible.

Options for the TE II include corrosion-resistant trough covers, storage cubicles and annex assemblies that connect the TE II to countertop work surfaces.

For detailed information about the TE II and other Sheldon products, call (601) 892-2731 or visit sheldonlabs.com.

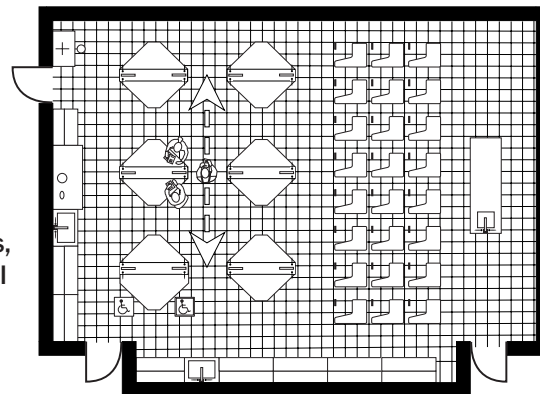

P.O. BOX 836 102 KIRK STREET
CRYSTAL SPRINGS, MS 39059
TEL 601 892 2731 FAX 601 892 4364
www.sheldonlabs.com